 **Swiss Diamond**®

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**INSTRUCTION MANUAL
ELECTRIC TOOTHBRUSH
SD-STB54805W**

ELECTRIC TOOTHBRUSH

Thank you for purchasing a product manufactured with Swiss Diamond standards. Our products are made in accordance with the highest demands of a quality, utility and design. We hope that you will continue to choose our products in the future.


Before using the product, make sure to read this manual carefully, it contains an important safety information and recommendations for use and correct maintenance.

Please save the instruction manual with a warranty card, cash receipt and if possible, with a gift box.

Introduction:

1. Description
2. Toothbrush head types
3. Important safeguards
4. Instruction for use
5. Features and properties
6. Charging of the battery
7. Cleaning and maintenance
8. Specification
9. Warranty
10. Certification

1. DESCRIPTION

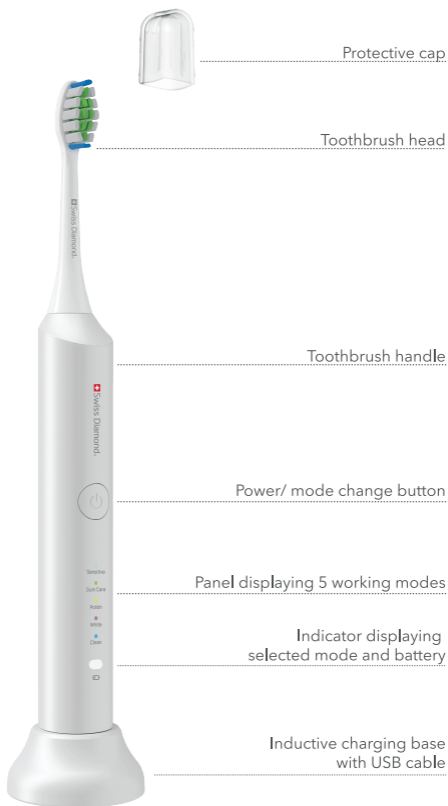2. TOOTHBRUSH HEAD TYPES

Replaceable toothbrush head for restoring the natural white color of the teeth. Effective at removing both the residue and discolorations.

Included.

Replaceable toothbrush head for tongue allows to clean the most hard-to-reach places surface of the tongue. Highly efficiency for cleaning teeth and tongue prevents accumulation of bacteria and guarantees fresh breath.

Included.

Replaceable toothbrush head for clean complex care of teeth and entire oral cavity. Full protection and fresh breath guaranteed every day.

Replaceable toothbrush head for gentle brushing of sensitive teeth and gums. Full comfort and maximum effectivity of brushing guaranteed.

CAUTION! It is recommended to replace the heads every 3 months with regular use.

3. IMPORTANT SAFEGUARDS

Read all instructions and important safeguards to avoid damage or injury when using this electric toothbrush:

1. The device should not be used in any way other than the one specified in the manual. The toothbrush is made only for cleaning teeth, gums and tongue.
2. Before plugging the device, make sure that the voltage shown on the device is the same as the one in the outlet.
3. If the device is damaged in any way or makes unusual noises after being powered on (in the area of head, main body or charger), do not use it. It is forbidden to try to repair it on your own. The device does not contain any elements made for the user's repairs.
4. The device should not be used by the people (including children) with physical, sensory or mental impairments and those with absence of knowledge or experience necessary to operate such devices, unless they are being supervised or instructed about usage of device by the person responsible for their safety.
5. The children should be instructed to not use the device as a toy.
6. When turning off the device using an adapter (not included), always hold the adapter itself, not the charger cord.
7. The inductive charger with USB. should be used and stored away from the heat sources, water and damp.
8. The manufacturer does not take responsibility for the damages caused by an improper use of the product.
9. Do not use the heads other than the ones recommended by the manufacturer.
10. Do not use the charger while in bath or under the shower.
11. Do not use this device in places where it can land or be thrown into the water. If the toothbrush has been dropped on

the floor, replace the head before next use, even if there is no visible damage.

12. Do not submerge the handle or the charger in water or any other liquid. After cleaning and before plugging the device, make sure that every component is dry.

13. The device is for the home use only. It is not made for professional use.

14. The device is for oral care only. Do not use for other purposes. Do not chew or gnaw on the heads.

15. If you are undergoing treatment for any oral cavity disease or had jaw/gum surgery within the last 2 months, make sure to contact your dentist before using the device.

16. Consult your dentist if your gums start bleeding or you feel pain after using the toothbrush.

17. If the user has a pacemaker or similar device, consult a dentist or product manufacturer before using the device.

18. Replace the heads every 3 months or sooner (if you notice the signs of wearing out).

19. Do not clean the heads in the dishwasher.

20. To prevent cracking, rinse the toothbrush head thoroughly with soap and water after each use, as peroxide, baking soda or other bicarbonates are often used in whitening toothpastes.

21. Only original spare parts may be used to repair the device.

22. Do not use the device after a fall if there are visible signs of damage. Before use, the device must be checked by a qualified person.

23. Always unplug your electric toothbrush when not in use.

24. Keep the power cord away from hot or sharp surfaces.

Avoid damage to the insulation of the power cord.

25. After cleaning the device and before connecting it to a power source, please make sure that the charging base is completely dry.

26. Use only original heads of this brand.

27. The use of accessories that are not included in the set of delivery will void your warranty.

28. Do not place metal objects on the charging base during charging, to avoid danger and damage to the device.

CAUTION! The device has a waterproof handle and is protected from moisture, but it is not recommended to use it under water!

4. INSTRUCTIONS FOR USE

01 Assembly

Before using the device, mount the brush head on top of the handle. There should be a small gap between the handle and the head. Avoid applying too much force while attaching the head to the handle. Do not twist the brush head, otherwise it will damage the pin and the motor.

02 Rinsing

Before every brushing, rinse the bristle with running water and apply the right amount of toothpaste on it.

03 Switching on the device

Press the On/Off button to switch the device on. With one press of the button, the toothbrush will start working.

04 Mode selection

The toothbrush offers five cleaning modes, allowing you to adjust it to your needs. Pressing the power button once (when the device is turned on) changes the operation mode in the following order:


- **Sensitive**
31000 rpm: Gentle brushing of sensitive teeth and surfaces.
- **Gum care**
35000 rpm: Massages the gums and improves their blood circulation
- **Polish**
41000 rpm: Tooth polishing and brightening
- **White**
45000 rpm: Tooth whitening, removal of stains and discoloration
- **Clean** 48000 rpm: Removes tartar and bacteria

Advice: it is recommended to choose the sensitive mode during your first usage.

05 Cleaning

The bristle should always be inclined at a 45 degree angle, in order to clean every spot within our mouth cavity.

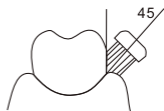
Advice: It is recommended to turn the device on after putting it in the mouth first.


06 Cleaning technique

Gently hold the handle of the toothbrush and move slowly in a circular movements, do not press hard, do not apply tension. During the cleaning, conditionally divide the mouth into 4 quarters for complete coverage of the oral cavity.

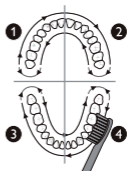07 Final rinsing

Wash your mouth out with water, then do the same with the toothbrush, removing bacteria and the remaining toothpaste.

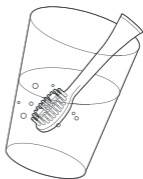5. FEATURES AND PROPERTIES

Waterproof level: IPX7

High waterproof level for comfortable and safe use.

One button

The device is operated by a single button. Pressing the button continuously after turning the toothbrush on allows us to change the cleaning mode.

Remembering the last mode chosen

After turning the device on again, it will be set to the previously chosen cleaning mode. It is necessary that the toothbrush operates in one selected mode for at least 5 seconds and then after turning it off, the device will remember the last brushing mode used.

Thirty-second interval

The toothbrush briefly stops working every 30 seconds, in order to let us know how much time has passed. This allows us to keep track of brushing every quarter of our oral cavity every 30 seconds. For this conditionally divide the oral cavity into 4 parts. Start cleaning each part in turn.

Two-minute timer

toothbrush can turn itself off automatically after 2 mins of usage. This function helps us to brush our teeth for two minutes, just like the dentists recommend.

Up to 240 minutes of work

After a single charge, the toothbrush can be used for 240 minutes (work in Sensitive mode).

6. CHARGING OF THE BATTERY

1. Before using the toothbrush for the first time, make sure that the device is fully charged.
2. Blinking of the battery indicator means that the battery power is running low.
3. Place the toothbrush on the inductive charging base to charge the battery. Connect the USB cable to an adapter (not included) and plug it into the power outlet. The charging start is indicated by flashing red.
4. A full charging cycle lasts from 8 to 16 hours allowing you to use the toothbrush for 60 days or less with regular use twice a day (depending on the selected mode).
5. When the battery is completely discharged, the toothbrush will stop working. After you plug the power adapter of the charging base with a discharged toothbrush into the mains, it may take 6 to 20 minutes before the toothbrush turns on and the indicator lights up. It is recommended to use a toothbrush with a fully charged battery. It will increase the toothbrush usage time.

WARNING! The device can only be used when it is charged! Mains usage (during charging) is impossible.

WARNING! When the battery is low, the charge indicator will not be as bright and the vibration intensity of the toothbrush will noticeably decrease, then you need to charge the battery.

7. CLEANING AND MAINTENANCE

1. Turn off the toothbrush and remove the replacement brush head after use. Carefully clean the top of the handle and pin.
2. After each use, rinse the toothbrush head and handle with warm running water. If necessary, you can use a soft cloth and soap to clean the surface.
3. Do not use organic solvents and abrasive materials for cleaning, as well as sharp objects that can damage the toothbrush coating.
4. Store your device in a cool and dry place.

Environment friendly disposal

Please remember to respect the local regulations: hand in the non-working electrical equipment to an appropriate waste disposal center.

8. SPECIFICATION

Product type:	Electric sonic toothbrush SWISS DIAMOND
Working modes:	Clean, white, polish, gum care, sensitive
Switch type:	A single button
Run time:	240 minutes (using sensitive mode)
Charging time:	8-16 hours
Battery:	Li-ion, 600 mAh, 3.7V
Battery indicator:	Flashing light
Noise level:	≤68 dB
Waterproof level:	IPX7
Power:	3 W
Voltage:	5V, 1A

9. WARRANTY

- 36 months of warranty is valid for this device.
- During the warranty period, the manufacturer is obliged to eliminate any manufacturing defects by repairing, replacing parts or replacing the entire product.
- The warranty comes into force only if the date of purchase is confirmed by the store's seal on the warranty card.
- This warranty is valid if the product has been used in accordance with the instruction manual, has not been repaired, disassembled or damaged as a result of improper usage.
- This warranty does not cover to body trim details, natural wear, consumables and accessories (heads, protective caps, etc.).

The full list of service centers you may find on the web-site: www.inhousewares.com

Life time is 5 years.

The month and year of manufacture are indicated on the bottom of the device and packaging.

The manufacturer has the right to make changes in the design, configuration, packaging, and also in the technical features of the product in order to improve products without further notice of these changes.

10. CERTIFICATION

This product meet the requirements of the technical regulations of the Customs Union, all required European and Russian safety and EMC standards.

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